



Peace of mind through accurate diagnosis

CHIRON® RIBA™ HCV 3.0 SIA

Hepatitis C Virus Encoded Antigen (Recombinant c33c and NS5 antigens; Synthetic 5-1-1, c100 and c22 peptides)

Globally recognized, CHIRON® RIBA™ HCV 3.0 SIA is intended to be used as a supplemental test for human serum or plasma specimens found to be repeatedly reactive in HCV antibody screening procedures. Accurate and reliable RIBA™ HCV 3.0 SIA test results will enable physicians and donor counselors to provide peace of mind to patients in determining their HCV status.

Assists in Donor Counseling

There is now significant evidence that an individual who is reactive for anti-HCV in a supplemental assay is likely to be infected with HCV.

— FDA Guidance for Industry, September 1998

RIBA™ HCV 3.0 SIA will prevent thousands of Unnecessary lookback notifications and subsequent testing and counseling of blood recipients based on inconclusive (indeterminate) HCV test results.

— FDA Talk Paper, February 12, 1999


Aids in Diagnosis

Supplemental anti-HCV testing confirms the presence of anti-HCV which indicates past or current infection...

— MMWR: 16 October 1998; 47(RR-19): 27.

CHIRON® RIBA™ HCV 3.0 SIA

offers superior performance

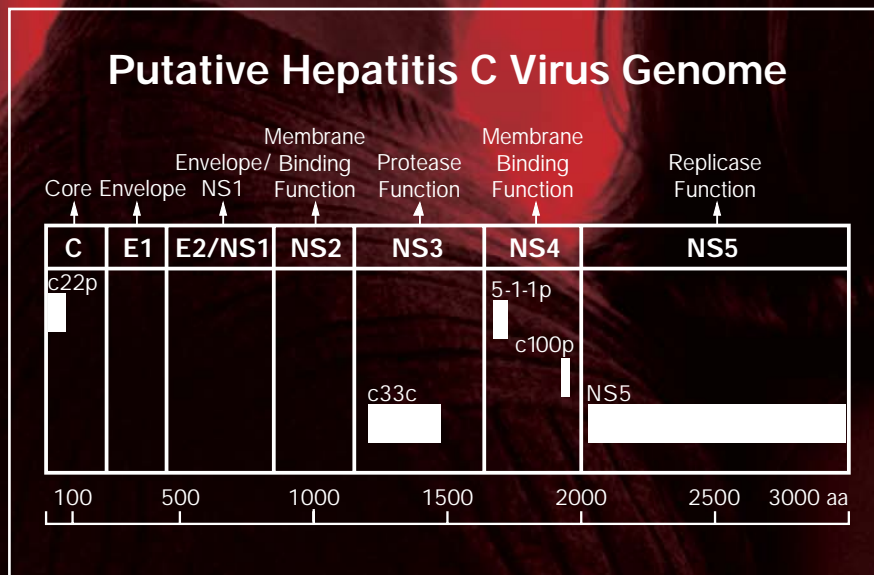


Improved Seroconversion Sensitivity

- In 52 out of 86 (60%) seroconversion panels, RIBA™ 3.0 SIA showed sensitivity or earlier positivity.

High Correlation with HCV RNA

- In a study of 220 HCV RNA positive individuals, 99.5% were reactive using RIBA™ 3.0 SIA.
- *Damen et al. Transfusion 1995; 35(9): 745-749.*

Identity of CHIRON® RIBA™ HCV 3.0 SIA Antigens

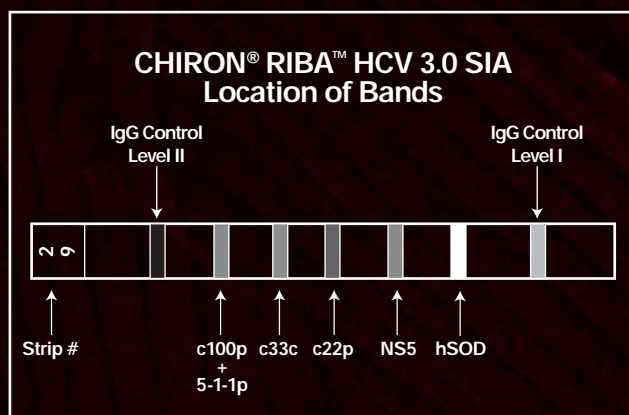

Genetically engineered HCV antigens and synthetic peptides improve performance

Ordering Information

| Product Description | Package Size | Product Code |
|---------------------------|--------------|--------------|
| CHIRON® RIBA™ HCV 3.0 SIA | 30 Tests | 930600 |

To place an order, please contact Customer Service at 1-800-828-6316 or order online at www.orthoclinical.com

Manufactured by

CHIRON

Chiron Corporation
4560 Horton Street
Emeryville CA 94608-2916

www.orthoclinical.com

Ortho® is a registered trademark of Ortho-Clinical Diagnostics.
CHIRON® RIBA™ HCV 3.0 SIA is a registered trademark of Chiron Corporation.

© Ortho-Clinical Diagnostics, Inc. 2002 DS1010 7/02

Distributed by

Ortho-Clinical Diagnostics

a Johnson & Johnson company

1001 US HWY 202
Raritan NJ 08869-0606